cAMP complete ELISA kit

Most sensitive (0.006 pmol/ml and 0.027 pmol/ml) and complete colorimetric kit for intra- and extracellular cAMP quantification.

This is a colorimetric competitive immunoassay kit for the quantitative determination of extracellular cAMP diluted in buffer or intracellular cAMP in cells or tissues lysed in 0.1M HCl. 0.1M HCl is used to stop endogenous phosphodiesterase activity and stabilize the released cyclic AMP. The 0.1M HCl treated samples are then analyzed directly in a microtiter plate without extraction, drying and reconstitution. Sensitivity is increased >10-fold by acetylation (reagents included). Absorbance is read at 405 nm. The large signal-to-background ratio offers superior sensitivity compared to competitors. This kit lets you easily assess adenylyl cyclase activation by GPCRs.

Cyclic AMP (cAMP) is one of the most important “second messengers” involved as a modulator of physiological processes. A number of hormones are known to activate cAMP through the action of the enzyme Adenylate cyclase which converts ATP to cAMP. cAMP has been shown to be involved in the cardiovascular and nervous systems, immune mechanisms, cell growth and differentiation, and general metabolism.

Citations: 47 View Online »

Ordering Information


Order Online »

ADI-901-163A	5x96 wells


Manuals, SDS & CofA

View Online »


- Versatile and convenient kit used to assay both intracellular and extracellular cAMP
- Highy sensitive measurement of cAMP with an optional acetylation protocol increasing sensitivity >10-fold (0.006 pmol/ml and 0.027 pmol/ml) for intra- and extracellular cAMP, respectively
- Simple, efficient and well-established sample handling protocols
- Widely cited in peer-reviewed publications
- Higher throughput format with results in 3 hours for up to 39 samples in duplicate
- Reliable and consistent lot-to-lot performance
- Save with bulk kit package!

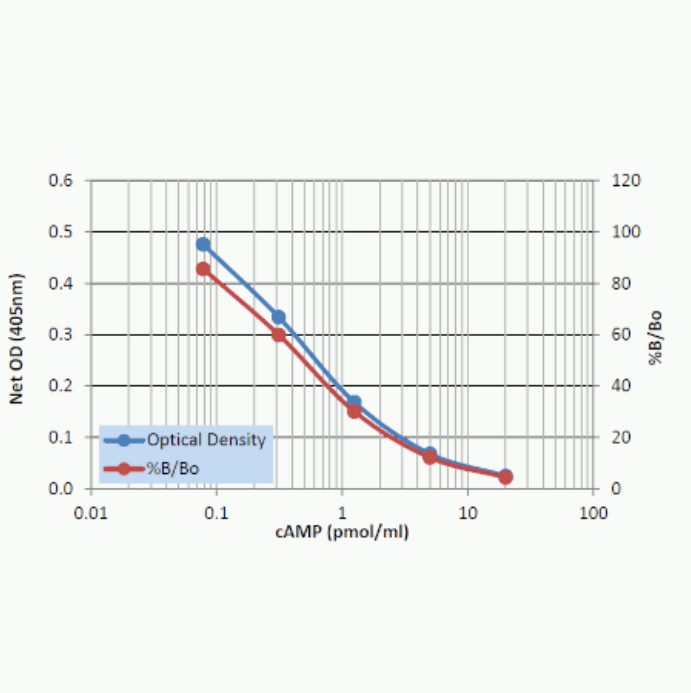


Typical standard curve used to calculate cAMP concentrations in 0.1M HCl for Acetylated assay format.

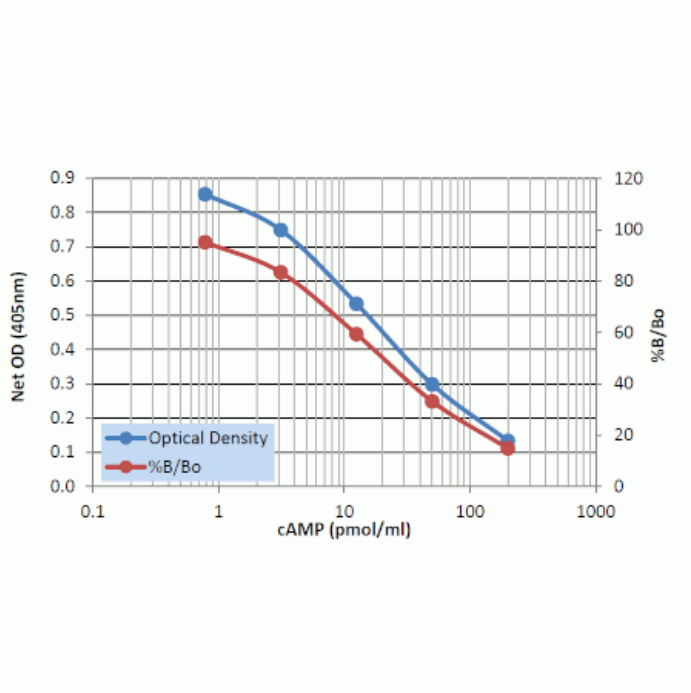

Typical standard curve used to calculate cAMP concentrations in 0.1M HCl for Non-acetylated assay format.

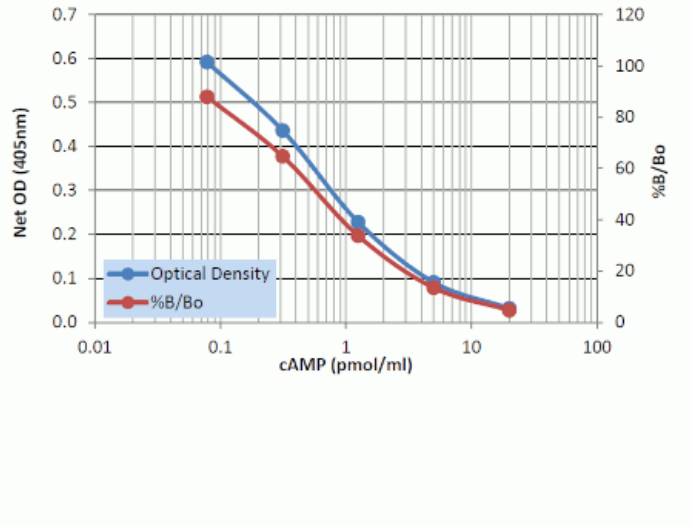

Typical standard curve used to calculate cAMP concentrations in Assay Buffer 2 for Acetylated assay format.

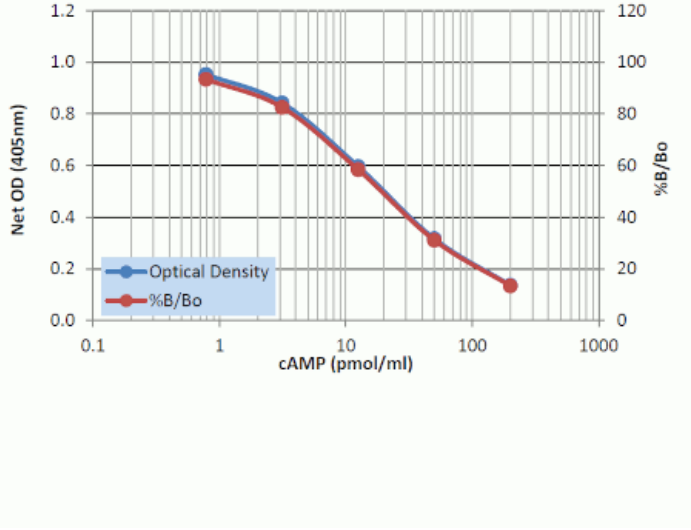

Typical standard curve used to calculate cAMP concentrations in Assay Buffer 2 for Non-acetylated assay format.

Handling & Storage

Use/Stability	Store all components at +4°C, except standard and conjugate at -20°C.
Shipping	Blue Ice

Regulatory Status

RUO - Research Use Only

Product Details

Alternative Name	cyclic AMP
Application	Colorimetric detection, ELISA
Application Notes	For the quantitative determination of cAMP in cell lysates, culture supernatants, serum, saliva, and tissue from any species.
Assay Time	3 hours
Compatibility	This product is compatible with the **Absorbance 96 Plate Reader**.
Contents	GxR Microtiter plate, Conjugate, Antibody, Assay buffer 2 concentrate, 0.1 M HCl, Neutralizing reagent, Wash buffer concentrate, Standard, pNpp Substrate, Stop solution, Acetic anhydride, Triethylamine
Crossreactivity	cyclic AMP (100%)AMP, ATP, cyclic GMP, GMP, GTP, cyclic UMP, CTP (<0.001%)
Sensitivity	Non-acetylated: 0.49 pmol/ml (assay buffer) / 1.18 pmol/ml (HCl) (range (non-acetylated) 0.78 – 200 pmol/ml)Acetlylated: 0.027 pmol/ml (assay buffer) / 0.006 pmol/ml (HCl) (range (acetylated) 0.078 – 20 pmol/ml)
Species Reactivity	Species independent
Wavelength	405 nm

ENZO LIFE SCIENCES, Inc.
Phone: 800.942.0430
info-usa@enzolifesciences.com

European Sales Office
ENZO LIFE SCIENCES (ELS) AG
Phone: +41 61 926 8989
info-eu@enzolifesciences.com

Belgium, The Netherlands & Luxembourg
Phone: +32 3 466 0420
info-be@enzolifesciences.com

France
Phone: +33 472 440 655
info-fr@enzolifesciences.com

Germany
Phone: +49 7621 5500 526
info-de@enzolifesciences.com

UK & Ireland
Phone (UK customers): 0845 601 1488
Phone: +44 1392 825900
info-uk@enzolifesciences.com

For local distributors and detailed product information visit us online: enzo.com All product names, logos, and brands are the property of their respective owners.